



**User manual
for
Roach formula**

Version 1, December 2023, in English

OBSOLETE

1. The Evidencio platform

The Evidencio platform facilitates the creation, use, validation and implementation of medical prediction models and clinical decision support tools. This user manual specifically relates to the Roach formula. The User Manual can also be referred to as the Instructions For Use (IFU).

Throughout this manual CE-marked content and the term medical device are used interchangeably.

2. Disclaimer

Evidencio provides information, models, calculators, equations, and algorithms (tools) intended for use by healthcare professionals. Some of these tools have been certified as CE-medical devices. For such CE-marked content the 'Official Legal Disclaimer for CE-marked content' applies. All other content and tools provided by Evidencio are explicitly only covered by the 'Official Legal Disclaimer for non CE-marked content' both are available here:

<https://www.evidencio.com/disclaimer>

3. Warnings



1. Warnings for CE-marked content

Calculations alone should never dictate patient care, and are no substitute for professional judgement. This tool is only to be used by physicians in a clinical setting, and is not for patient use.

Always read the intended use before using this tool.


Before reading the result, double check the filled in values to prevent errors.

Results that concern risk percentages, do not guarantee certain outcomes. When there is a risk present, do not expect an event to not occur at all, even if the risk is very small.

This model is only intended for use in settings where the usage and result of a model are never immediately needed.

4. Model landing page

The medical device model on the Evidencio platform is shown in Figure 1. The model landing page contains the following sections, that are indicated in Figure 1.

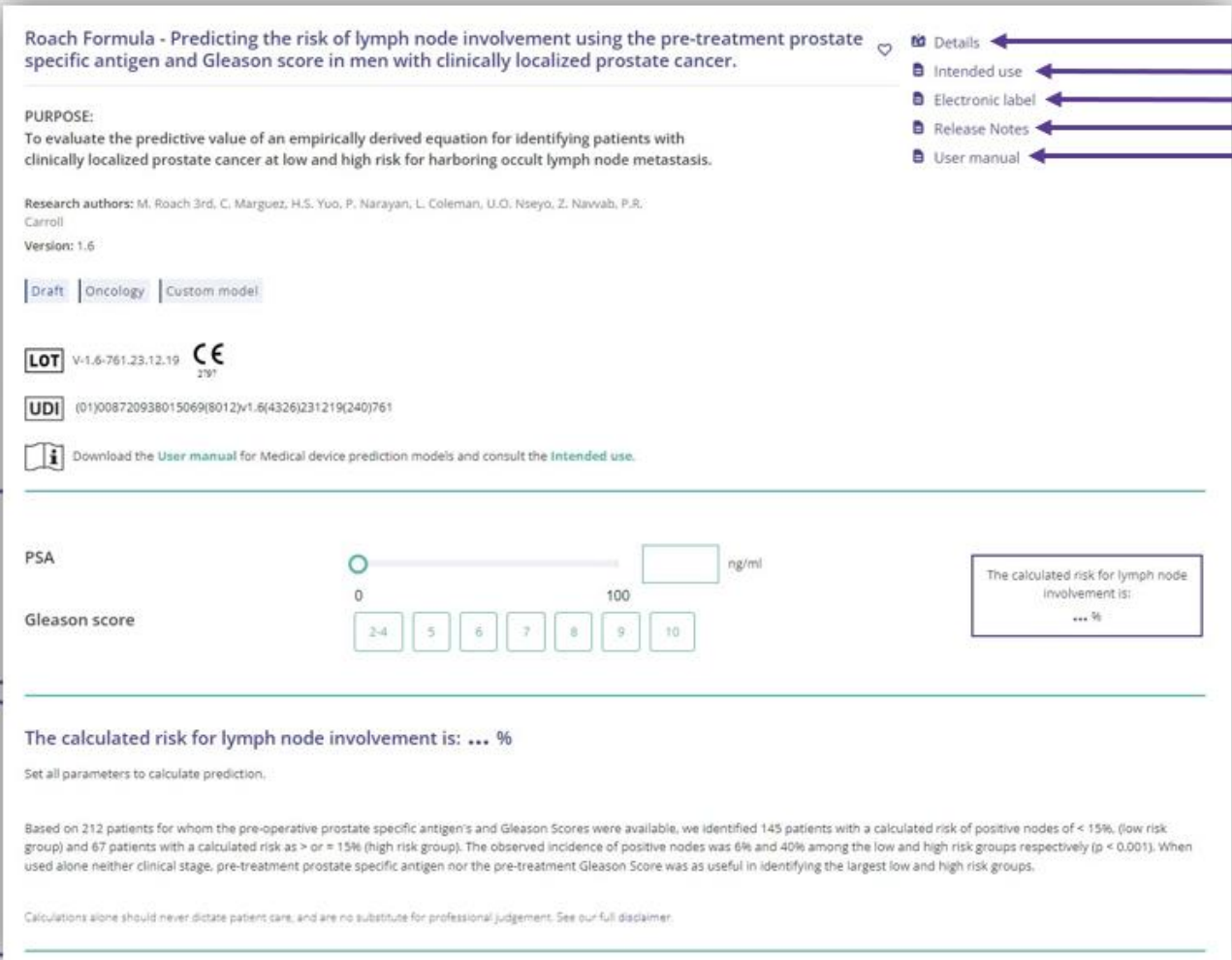

A. Model title
Roach Formula - Predicting the risk of lymph node involvement using the pre-treatment prostate specific antigen and Gleason score in men with clinically localized prostate cancer.

B. Model description
PURPOSE:
To evaluate the predictive value of an empirically derived equation for identifying patients with clinically localized prostate cancer at low and high risk for harboring occult lymph node metastasis.

C. Research authors
Research authors: M. Roach 3rd, C. Marguez, H.S. Yuo, P. Narayan, L. Coleman, U.O. Nseyo, Z. Navvab, P.R. Carroll
Version: 1.6

D. Model tags
Draft | Oncology | Custom model

E. LOT
V-1.6-761.23.12.19

F. UDI
(01)008720938015069(8012)v1.6(4326)231219(240)761

K. Download the User manual
Download the [User manual](#) for Medical device prediction models and consult the [Intended use](#).

L. Input fields
PSA: 0 to 100 ng/ml
Gleason score: 2-4, 5, 6, 7, 8, 9, 10

M. Calculated risk
The calculated risk for lymph node involvement is: ... %

Set all parameters to calculate prediction.

Based on 212 patients for whom the pre-operative prostate specific antigen's and Gleason Scores were available, we identified 145 patients with a calculated risk of positive nodes of < 15% (low risk group) and 67 patients with a calculated risk as > or = 15% (high risk group). The observed incidence of positive nodes was 6% and 40% among the low and high risk groups respectively (p < 0.001). When used alone neither clinical stage, pre-treatment prostate specific antigen nor the pre-treatment Gleason Score was as useful in identifying the largest low and high risk groups.

Calculations alone should never dictate patient care, and are no substitute for professional judgement. See our full disclaimer.

G. Details
H. Intended use
I. Electronic label
J. Release Notes
K. User manual

Figure 1. An example of a model landing page.

A. Model title

This is the title and name of the model.

B. Model description

This is a short description of the model.

C. Research authors

These are the research authors of the paper that originally published the model.

D. Model tags

These are the tags that are assigned to the model. Evidencio has the following status tags: "Draft", "Public", "Private", "Under review". Evidencio has the following model type tags: "Composite model", "Sequential model", "API model". Evidencio has the following calculation method tags: "Linear model", "Logistic regression", "Cox regression", "RScript" and "Custom model". Next to this, there are tags that indicate the specialty e.g. "Cardiology".

E. LOT number

The LOT number indicated the model version, the model identifier, and the model publication date. Publication date is indicated as YY.MM.DD.

Additionally, the CE mark is displayed next to the LOT number. This way, medical devices can be easily recognized.

F. UDI number

The UDI number is an international tool that helps users identify and find information on products. UDI stands for Unique Device Identifier. Evidencio's UDIs have the following format:

(01)UDI-DI number(8012)versionnumber(4326)releasedate(240)identificationnumber

The UDI-DI number is a unique numeric code. For each medical device of Evidencio, a unique UDI-DI is ascribed. This UDI-DI is used as an "access key" for information stored in a unique device identification database (UDID). Information on Evidencio's medical devices can be found by searching for the UDI-DI number in the following data base:

<https://gepir.gs1.org/index.php/search-by-gtin>




G. Details button

On the top right of the model page, several clickable buttons are displayed that show a pop-up when clicked. The first button opens a pop-up concerning additional information about the model. This pop-up has three sections: Details, Study characteristics and Supporting publications & related files.

Details

The first part of the additional information concerns the details of the model as shown in Figure 2.

Details

Model author	T. A. Hueting	Status	Draft
Model ID	761	Share	  
Version	1.6		
Revision date	2023-12-19		
Specialty	Oncology , Urology		
Model type	Custom model (Conditional)		
MeSH terms	<ul style="list-style-type: none"> Prostate Cancer Lymphadenectomy Gleason Score Prostate Specific Antigen 		

Condition

Formula

$$\left(\frac{2}{3}\right) \cdot \text{PSA} + (\text{Gleason score} - 6) \cdot 10 < 0$$

$$\left(\frac{2}{3}\right) \cdot \text{PSA} + (\text{Gleason score} - 6) \cdot 10 \geq 0$$

$$\left(\frac{2}{3}\right) \cdot \text{PSA} + (\text{Gleason score} - 6) \cdot 10$$

Figure 2. The model details.

Study characteristics

Below the 'Details section' the section labeled 'Study characteristics' provides information on the characteristics of the patient data used to derive and validate the model. Additional information is provided on the methods used to develop and/or validate the model.

An important part of the Study characteristics is the information on Supporting publications and related files. These sections can be found at the bottom of the Details-pop-up as shown in Figure 3.

Supporting Publications

Title or description

Predicting the risk of lymph node involvement using the pre-treatment prostate specific antigen and Gleason score in men with clinically localized prostate cancer.

Tags

Paper

Model coefficients

Related files

No related files available

Figure 3. An example of Supporting publications & related files.

H. Intended use button

Intended medical use

The device is intended to be used by healthcare professionals to estimate the probability of lymph node involvement in patients with clinically located prostate cancer. The device combines prostate specific antigen (PSA) and the Gleason Score (GS) to predict lymph node involvement.

The device is intended to be used for patients with clinically located prostate cancer, by healthcare professionals and qualified medical specialists in a clinical setting. The device is not intended for use by patients on their own.

The Roach formula is not intended to replace clinical decision-making, it can only inform the healthcare professional, and only provides a cumulative percentage to estimate the risk of lymph node involvement. No direct instructions for further diagnostics, treatment, or otherwise, are given.

Clinical Benefit

The Roach formula is intended to assist patients with relevant and specified clinical outcome parameters. Concretely, this is achieved by estimating a risk in order to support clinical decision making aimed at patients with clinically localized prostate cancer in order to support clinical decision making regarding pelvic lymph node treatment. Correct functioning of the Roach formula can result in these clinical benefits:

- The Roach formula can assist in risk stratification for patients.
- Risk stratification can reduce the burden of (invasive and intensive) medical procedures such as tests on patients with low risks, reducing, shortening or avoiding stays in hospitals or other care facilities.
- Risk stratification can reduce the unnecessary consumption of (scarce) medical resources, decreasing costs and increasing their availability for high risk patients.

Intended target population and exclusion

The Roach formula is intended to be used only for a specific group of patients, corresponding to the below indications and contra-indications. Additionally, it should be noted that patient data for derivation was gathered exclusively in patients that underwent a radical prostatectomy for prostate cancer.

Inclusion Criteria

The Roach formula is intended for:

- Patients with clinically localized prostate cancer, eligible for radical prostatectomy

Exclusion Criteria

The Roach formula should not be used for patients that do not fulfil the inclusion criteria. It is not intended for:

- Patients with incomplete pathologic or biopsy data

User profile

Since Lymph Node Involvement in Prostate Cancer is regarded as a 'critical healthcare situation or condition', the use of the SaMD is intended for healthcare professionals and qualified medical specialists in a clinical setting. Healthcare professionals do not require additional training prior to the use of the medical device. The SaMD should not be used by patients.

Intended use environment

The SaMD can be used as made available on the Evidencio platform in any actively supported web-browser on personal computers, mobile devices, or tablet PCs, and on the mobile app provided by Evidencio. Furthermore, the SaMD can be used through the Evidencio iFrame representation of the SaMD, as an embedded view, provided that the specific Evidencio guidelines for iFrame implementations of this SaMD are adhered to. In addition, automatic calculation for the

model can be used through Evidencio's API. The model is only intended for use in settings where the usage and result of a model are never immediately needed.

Functioning, physical principle

The SaMD's underlying mathematical formula is a custom model. The acquisition and processing of the data, the analyses to assemble the relevant criteria for the SaMD as well as the setup and refinement of the Roach formula are described in the original studies from Roach M 3rd et al. Entering the details of an individual in the Evidencio SaMD initiates the calculation of the estimated risk of lymph node involvement.

I. Electronic label button

The electronic label button opens a pop-up with the location and address of Evidencio, the LOT number, the UDI number, the CE-mark, the medical device logo and a download link for the declaration of conformity of the medical device. The example of the electronic label is shown in Figure 4.

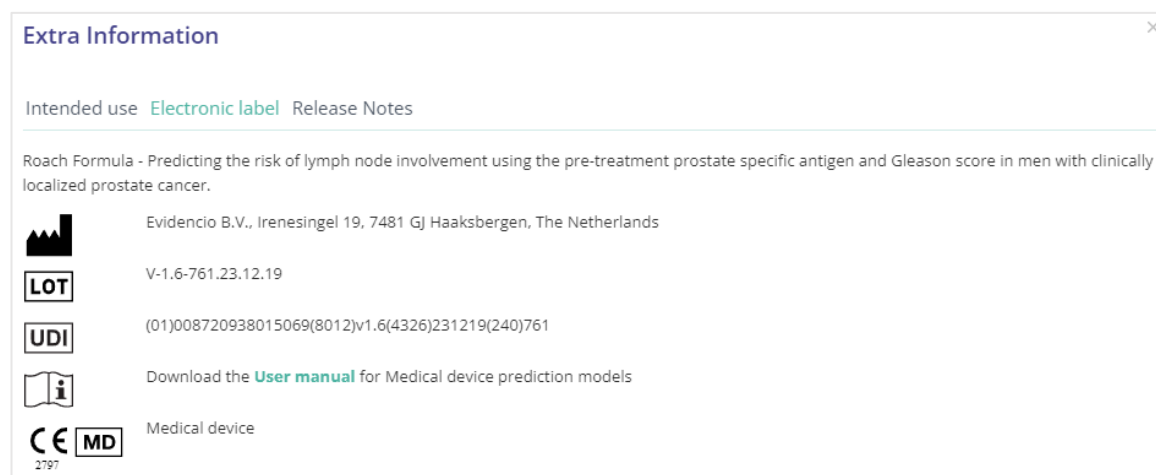


Figure 4. Example of the electronic label

J. Release notes

The 'Release Notes' button opens a pop-up with the latest release notes of the model. Here you can find what has changed over the last versions of the model. Additionally, if there are any known residual anomalies the user should be aware of, they are listed here.

K. User Manual

This user manual can be found in three places: 1) under the short description, 2) on the right of the model page, and 3) in the electronic label. Additionally, all versions of the user manual can be found in the general page for all user manuals for medical devices. The page can be found under the 'About' drop-down menu button as shown in Figure 5. The user manual page is shown in Figure 6.

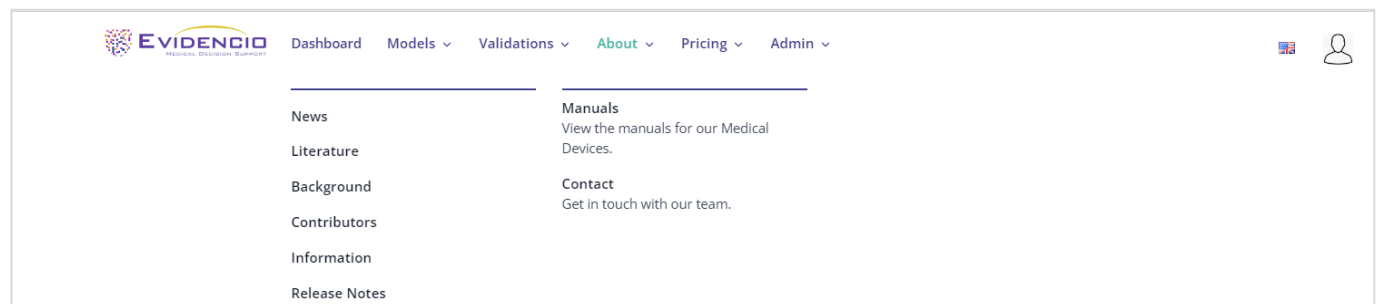


Figure 5. The drop-down menu where the user manual page can be found.

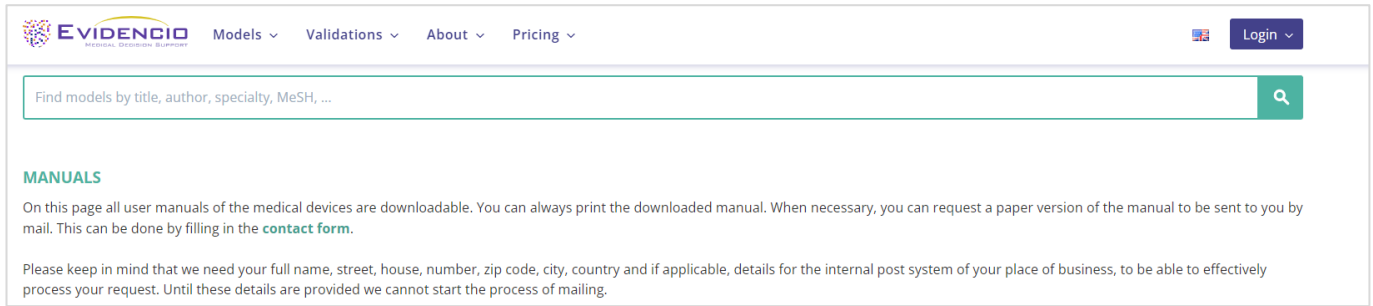


Figure 6. The user manual page for all user manuals.

You (The user) can always print this downloaded manual. When necessary, you can request a paper version of the manual to be sent to you by mail. Evidencio's contact details are listed in Chapter 6 of this user manual.

L. Input section

The Evidencio platform allows two separate input variables; categorical, and continuous variables.

Categorical variables

In the example shown in Figures 7 and 8, the **Gleason score** variable concerns a categorical variable. The patient status can be entered by clicking on either button. The selected button changes to green, as seen in Figure 8.

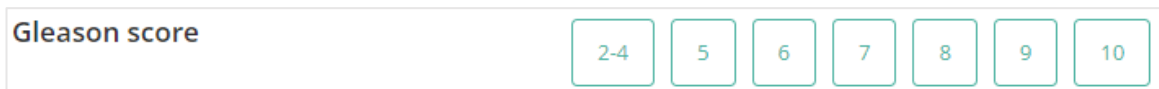


Figure 7. The variable for Gleason score, where no button has been clicked, and thus no input has been provided by the user.

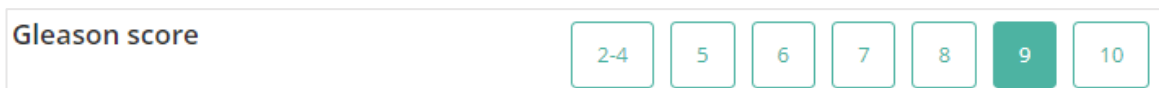


Figure 8. The variable for Gleason score, where the "9" button has been clicked.

Continuous variables

In the example shown in Figure 9, the **PSA** variable, exemplifies a continuous variable. The plausible ranges for the variables are used for the model.

The details for a patient can be entered by sliding the button to the correct value, or by entering the correct value in the box on the right-hand side (i.e., where the 20 is entered for **PSA**).

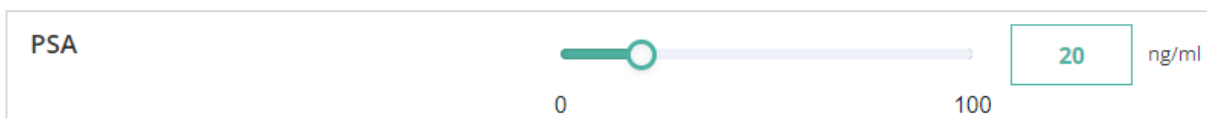


Figure 9. The variable for PSA, where "20" has been entered

M. Result section

At the bottom of the page, the results of the model are shown.

Result calculation

When all variables are filled in, a result will be calculated. No risk is displayed until all variables are filled in. The result section indicates "Set all parameters to calculate prediction."

Result interpretation

In the result interpretation, a risk stratification is given based on the risk score. An example of the information is shown In Figure 10.

The calculated risk for lymph node involvement is: **43 %**

Based on 212 patients for whom the pre-operative prostate specific antigen's and Gleason Scores were available, we identified 145 patients with a calculated risk of positive nodes of < 15% (low risk group) and 67 patients with a calculated risk as > or = 15% (high risk group). The observed incidence of positive nodes was 6% and 40% among the low and high risk groups respectively ($p < 0.001$). When used alone neither clinical stage, pre-treatment prostate specific antigen nor the pre-treatment Gleason Score was as useful in identifying the largest low and high risk groups.

Calculations alone should never dictate patient care, and are no substitute for professional judgement. See our full disclaimer.

Figure 10. The result information

Relevant information for correct use of the model

At the bottom of the page, there is a link to Evidencio's terms and conditions of use, the privacy policy, and the disclaimer.

5. Use of Medical devices

In general, and unless explicitly stated otherwise, CE-marked tools on Evidencio are only to be used by physicians in a clinical setting, and are not for patient use.

To use the tool, Evidencio requires a stable internet connection and runs on the following devices:

- Personal computers or laptops using the following browsers:
 - Safari (the latest three versions)
 - Chrome (the latest three versions)
 - Firefox (the latest three versions)
 - Edge (the latest three versions)
- Tablets or smartphones running on the next operating systems:
 - IOS (the latest three versions)
 - Android (the latest three versions)

The medical device cannot be used in combination with Internet Explorer. The personal computers, laptops, tablets or smartphones used should at least be able to have an internet connection and use the browsers mentioned above. The minimal screen resolution should be 800x600.

Furthermore, the model may be used through the Evidencio iFrame representation of the calculator, as an embedded view, provided that the specific Evidencio guidelines for iFrame implementations of that model are adhered to.

The Evidencio SaMD models can be used with any browser settings that don't distort the regular display of websites, with a 50% to 500% zoom rate, and at a display resolution starting from 800x600. However, factory recommended browser settings, 100% zoom rate and regular display resolution are recommended.

This model is only intended for use in settings where the usage and result of a model are never immediately needed.

6. Manufacturer details

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the country in which you, the reader, are established. A competent authority is the institute that governs all issues related to medical devices in a country.

Contact details of your competent authority can be found here: <https://www.ema.europa.eu/en/partners-networks/eu-partners/eu-member-states/national-competent-authorities-human>

Please contact Evidencio when you suspect any malfunction or changes in the performance of a medical device. Do not use the device, until Evidencio replies to your message that it is safe to start using it again.

Contact details of Evidencio:



Evidencio B.V., Irenesingel 19, 7481 GJ Haaksbergen, The Netherlands
www.evidencio.com
tel: +31 53 85195 08
e-mail: info@evidencio.com